

Wrist Blood
Pressure Monitor

UB-533



Wrist Blood
Pressure Monitor

UB-525



Wrist Blood Pressure Monitor UB-533/UB-525

UB-533

All-in-one Wrist Monitor



The UB-533 is designed for 2 users and focuses on greater ease of use with professional quality. It is equipped with functions like CPG, cuff fit and a movement error indicators to support users in accurate measurement.

The %IHB and WHO blood pressure classification indicators provide users with detailed monitoring.

UB-525

Simple Operation
with Accuracy






The UB-525 is equipped with all functions essential for blood pressure monitoring. Its sleek design comfortably fits the wrist and provides accurate measurement.

About "CPG" Function(UB-533)



CPG (Correct Position Guidance) indicates the correct arm position with both the color and position of the light.

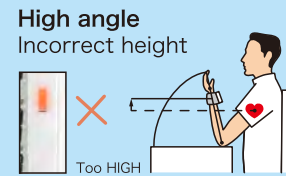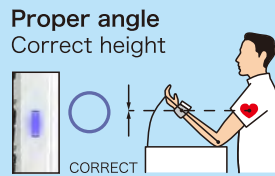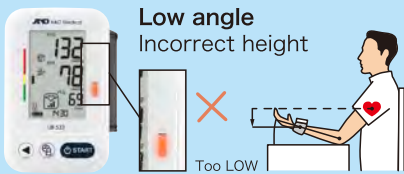

Specifications

	UB-533	UB-525
Measurement method	Oscillometric measurement	
Display type	Digital LCD, Pressure/Pulse displayed simultaneously	
Measurement range	Pressure : 0 - 299 mmHg Pulse : 40-180 beats / minute	
Measurement accuracy	Pressure: ± 3 mmHg Pulse: $\pm 5\%$	
Power supply	2 x 1.5V alkaline batteries (LR03 or AAA)	
Number of measurements	Approx. 250 measurements, when AAA alkaline batteries are used, with pressure value of 170 mmHg at room temperature of 23°C.	
Wrist circumference	13.5 - 21.5 cm	
Memory storage	60 x 2 users	60
Memory data average	Maximum 60 data average of systolic, diastolic, and pulse values	
Operating conditions	+10°C to +40°C / 15%RH to 85 %RH / 800 hPa to 1060 hPa	
Transport/Storage conditions	-20°C to +60°C / 10%RH to 95%RH	
Dimensions	Approx. 56 [W] x 88 [H] x 21.5 [D] mm	
Weight	Approx. 105g, excluding the batteries	



Correct Position Guidance



Cuff Fit Error Indicator



Movement Error Indicator



%IHB Indicator to show frequency of Irregular Heartbeat



WHO Classification Indicator



Irregular Heartbeat Indicator



2 user memory (60 x 2)



60-memory



Clock & Calendar



Average reading of last 60 measurements

Standard accessories

- Blood pressure monitor
- Instruction manual
- Batteries for temporary use x 2
- Carrying case(UB-533)

AND
A&D Medical

A&D Company, Limited

3-23-14 Higashi-Ikebukuro, Toshima-ku, Tokyo 170-0013 Japan
Telephone:[81](3)5391-6132 Fax:[81](3)5391-6148
1-243, Asahi, Kitamoto-shi, Saitama, 364-8585 Japan
Telephone:[81](0485)93-1111 Fax:[81](0485)93-1119
<http://www.aandd.jp>

A&D INSTRUMENTS LIMITED

Unit 24/26 Blacklands Way Abingdon Business Park,
Abingdon, Oxon OX14 1DY UNITED KINGDOM
Telephone:[44](1235) 550420 Fax:[44](1235) 550485
(Authorized Representative Established in the European Community)